



# TWO IN ONE!

VERSATILITY IS A BREEZE WITH THE P100

*The BREEZY P100 is a folding powered wheelchair suitable for indoor and outdoor use. It offers quick assembly and disassembly, and is easily converted to a manual chair with the simple addition of self propelling 24" wheels included in the package.*

## FROM POWER.....



### POWERED FUNCTION

Simply attach the portable motor unit under the base of the chair and you're ready to zoom around town.....

## .....TO MANUAL



### MANUAL FUNCTION

....OR remove the powered unit and simply attach the 24" Quick release PU Wheel provided.



# IN TWO EASY STEPS!



- **HYBRID** – both a power & manual chair in one
- **AUSTRALIAN STANDARDS** tested
- **ELECTRONICS** – Easy to use VR2 electronics
- **LITHIUM ION BATTERIES** for longer life expectancy
- **COLOURS** available in Red
- **VERSATILE** for indoor & outdoor use
- **FOLDING** convertable power wheelchair
- **BACK POUCH** for convenient storage
- **ADJUSTABLE** seat height

## FEATURES



**COMPACT & FOLDABLE**



**S/A DETACHABLE FOOTREST**



**FLIP BACK ARMRESTS**




**HALF FOLDING BACKREST**



**USER FRIENDLY ELECTRONICS**

## DIMENSIONS

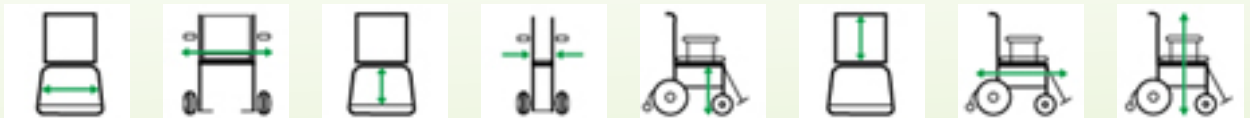

CODE	Seat Width	Frame Width	Seat Depth	Folding Width	Armrest Height	Back Height	Frame Depth	Adj Seat Height
BRP100-16	400mm (16")	520mm (20")	398mm (16")	300mm (12")	210mm (8")	435mm (17")	1040mm (41")	880 – 930mm (34.5" - 36")
BRP100-18	450mm (18")	570mm (22")	398mm (16")	300mm (12")	210mm (8")	435mm (17")	1040mm (41")	880 – 930mm (34.5" - 36")

## SPECIFICATIONS

<b>Front wheels</b>	200*50 PU Tyre	<b>Max. Groove Crossing Ability</b>	1000mm
<b>Rear wheels</b>	2.50 – 6 pneumatic tyre / Max. 50PSI (3.5 BAR)	<b>Max. Speed Forward</b>	7.8km / hr
<b>Seat Width</b>	450mm (18"), 400mm (16")	<b>Max. Speed in Low Speed Mode</b>	2km / hr
<b>Motor Output</b>	200W*DC24V*2PCS	<b>Max Weight Capacity</b>	120kg
<b>Brakes</b>	Magnetic Brake	<b>Weight (without batteries)</b>	33.4kg
<b>Battery</b>	20AH * 24V * 1PCS (Lithium Battery)	<b>Front Frame Weight</b>	15.4kg
<b>Range</b>	25km (on the flat road, 120kg load)	<b>Rear Frame Weight</b>	18kg
<b>Slope Grade Durability</b>	10° max.	<b>Battery Weight</b>	5.4kg
<b>Min Turning Radius</b>	Approx. 870mm	<b>Anti-tip Wheel</b>	Yes (on power unit only)
<b>Max. Kerb Climbing Ability</b>	40mm	<b>Colour</b>	Red