




The New Standard in Homecare



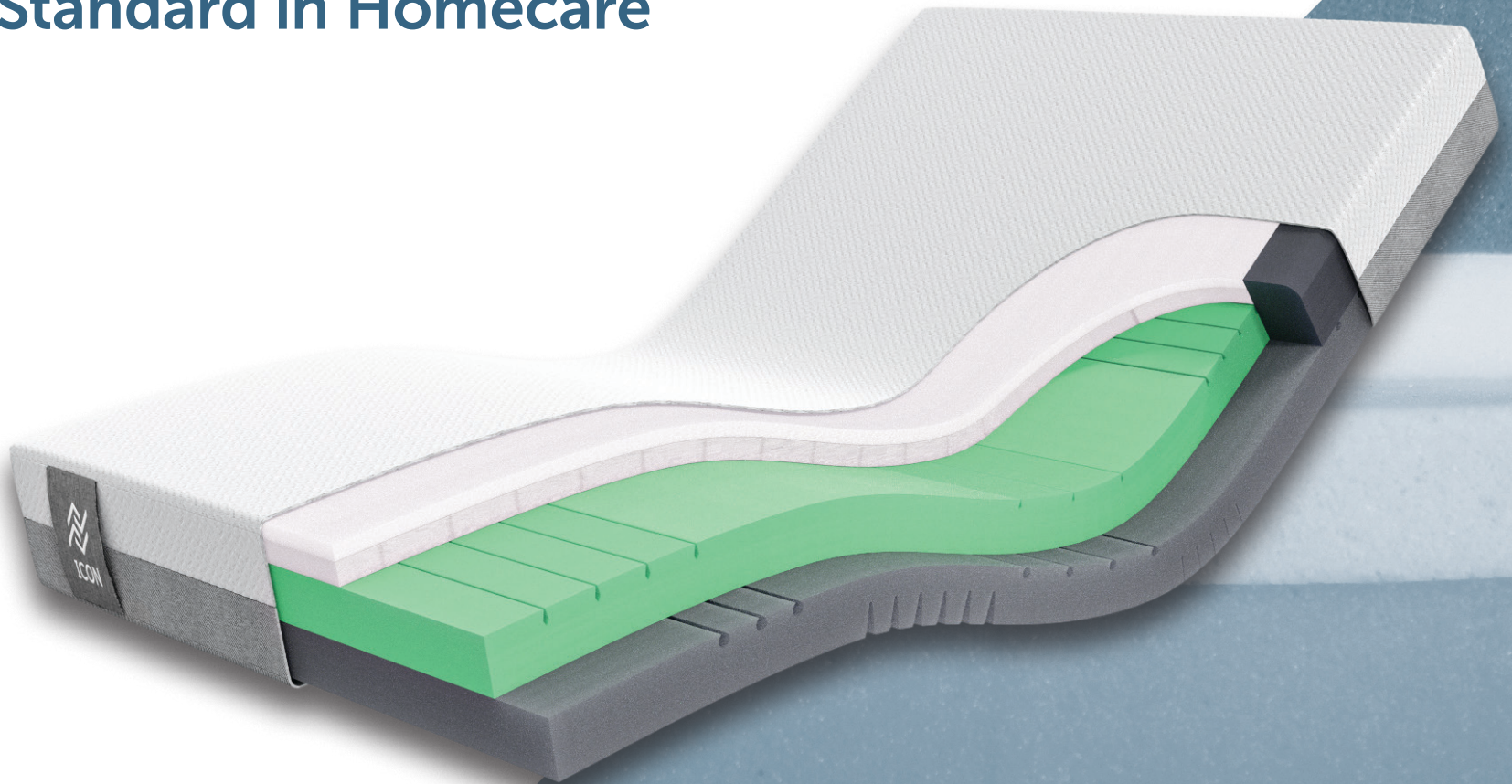

Designed for Safety, Engineered for Comfort

Compatibility with Adjustable Beds

- Compatible with a wide range of adjustable bed frames
- Custom thickness (15-20 cm) with consistent comfort

Breathable & Comfortable for Homecare

- Advanced airflow for temperature and moisture control
- Medical-grade support with a home-like feel

Easier Repositioning & Increased Independence

- Low-resistance design for easy movement
- Soft, Medium and Firm options for support needs

Safe & Supportive Transfers

- Contoured edges support easier entry and exit
- Stable surface for confident transfers



Risk Rating

	1	2	3	4	5
Shear & Friction	<div><div></div></div>				
Repositioning	1	2	3	4	5
	<div><div></div></div>				
Infection Control	1	2	3	4	5
	<div><div></div></div>				
Firmness	1	2	3	4	5
	<div><div></div></div>				
Moisture & Micro-climate	1	2	3	4	5
	<div><div></div></div>				
Performance & Longevity	1	2	3	4	5
	<div><div></div></div>				
Mobility & Function	1	2	3	4	5
	<div><div></div></div>				
Pressure Care Rating	Low - Medium Risk				
Warranty	5-year foam core 2-year Luxury Quilted Homecare cover				
Therapeutic Weight Range	Soft (150 kg) Medium (180 kg) Firm (300 kg)				

**Minimum and maximum weight capacities are dependent on client anatomy, bed profiling functions, and positioning.*



Product Features

Breathability & Microclimate Control

- Airoform Air Flow Interface enhances ventilation
- Humidity-resistant foam formulation for temperature regulation

Perimeter Support & Stability


- High-support perimeter for stability
- Supportive edges for transfers and repositioning

Superior Repositioning Support

- Designed for users with reduced mobility
- More stable and supportive than traditional reactive foam surfaces

Outstanding Comfort

- Dual-layer Airoform foam for effective pressure redistribution
- Available in multiple comfort options for individual needs

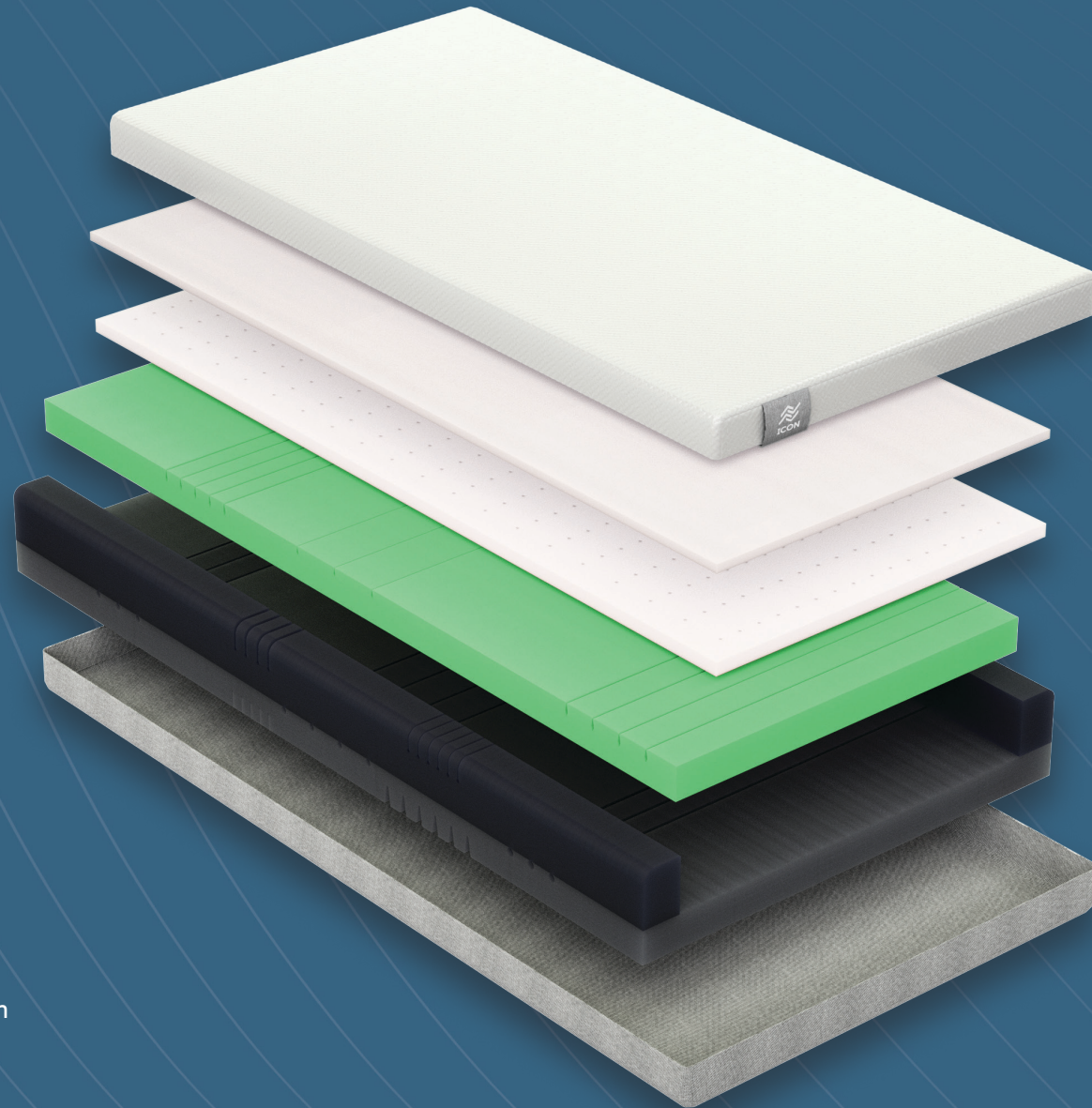

Foundational Support Foam

- Premium-grade firm foam enhances safety during transfers
- Maintains support and comfort in varying user positions

Advanced Medical Foam Technology

- Four-layer progressive support structure
- Engineered for durability and long-term performance

Adaptive Contouring for Adjustable Beds

- Unique profile hinge system for bed articulation
- Channelled profiling ensures a responsive and supportive feel

Hygiene & Longevity

- Cover options for different care environments
- Built with high-quality Advantiflex and side wall reinforcements

Cover Options

A range of high-performance covers tailored to suit different care environments and infection control needs.

Luxury Quilted Homecare

- 4-way stretch, anti-shear fabric
- High breathability, ideal for cooling
- Homely quilted finish with medical-grade performance
- Supports effective pressure care
- Removable and replaceable



Advantiflex Classic

- 4-way stretch, anti-shear surface
- Waterproof, wipe-clean finish
- High-frequency sealed and welded seams
- Durable PVC base for added protection
- Compatible with disinfectants up to 1,000 ppm




Premiflex Ultra

- 4-way stretch, anti-shear material
- Waterproof, breathable, and wipe-clean
- Enhances foam longevity and airflow
- Compatible with disinfectants up to 10,000 ppm



Standards of Safety & Care

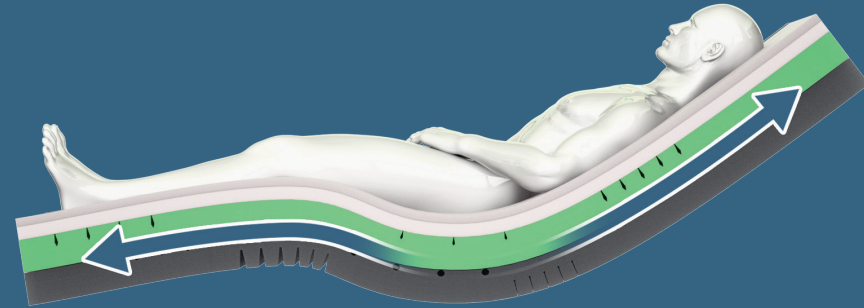

Breathability & Microclimate Control

- Profiling and channels enhance airflow for temperature regulation
- Open-cell foam structure promotes breathability
- Advanced cover design supports microclimate management
- Responsive foam top layer balances comfort and ventilation




Safer Transfers & Stability

- Radius side support promotes stability, allowing users to reposition with confidence
- Thick, high-density sidewalls enhance edge support for safer transfers
- Firm, supportive perimeter enhances stability.


Effortless Repositioning

- High-grade comfort layers make repositioning easier
- Reduced memory foam content for improved movement and clinical support
- Designed to differentiate from entry-level medical mattresses that restrict movement

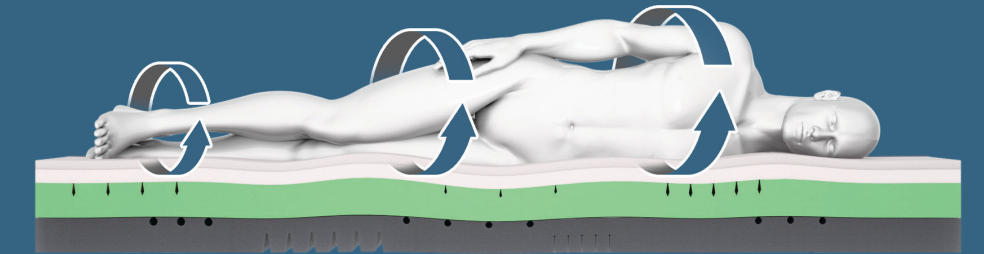

Universal & Adaptive Fit

- Adapts to movement, maximising adjustable bed benefits
- Knee and sacral profiling enhance flexibility
- Works with all comfort options and thicknesses
- Prioritises comfort without compromising safety

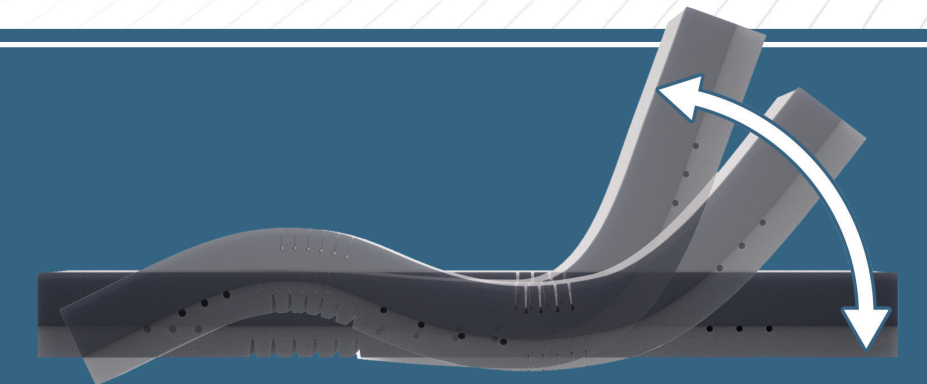


Comfort Options

The Icon mattress is available in Soft, Medium and Firm comfort settings, recognising that comfort is deeply personal and essential to effective clinical care. These options allow clinicians and caregivers to select a configuration that aligns with each client's unique needs, without compromising on the mattress's proven pressure care performance or core design principles.


ICON Soft

Ideal for clients who require gentle immersion and greater envelopment, helping reduce pressure peaks and enhance comfort for those with sensitivity or frailty.

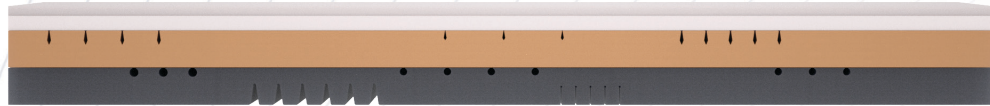


ICON Medium

A balanced option that suits a broad range of clients, offering both responsive support and comfortable contouring.

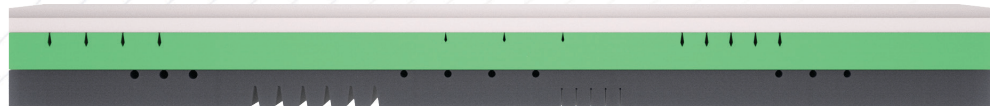

ICON Firm

Designed for clients who need added stability or prefer a more supportive surface, promoting easier repositioning and transfer support.

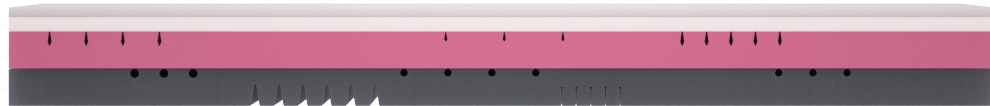


**Extra Soft and Extra Firm options are also available*

Height Options


The Icon mattress is offered in two height profiles, 15 cm and 20 cm, to provide flexibility in both clinical functionality and user preference. Each option has been designed to support specific care needs while maintaining the Icon's high standards in pressure care and durability.

15 cm Height

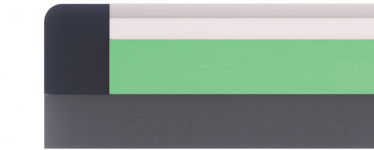

- Lower height option for adjustable bed configurations.
- Low-profile option for homecare environments
- Designed to align with key bed articulation points.
- Enables easier transfers for both clients and carers

20 cm Height

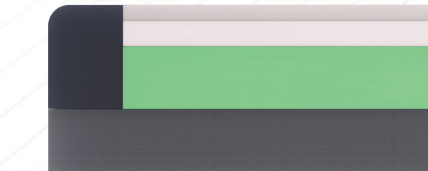

- Provides a more traditional mattress appearance and user experience
- Enhanced immersion and envelopment for improved pressure redistribution
- Available in Soft, Medium and Firm comfort options
- Increases long-term durability, ideal for high-use environments
- A versatile option for residential aged care and long-stay clinical settings

Custom Sizing Solutions

Specialist Mattress Manufacturers

Size	Dimensions (mm)	Soft		Medium		Firm	
		150 mm	200 mm	150 mm	200 mm	150 mm	200 mm
Single	1980x880	ICONS15S	ICONS20S	ICONS15M	ICONS20M	ICONS15F	ICONS20F
Single, Long	2030x880	ICONSL15S	ICONSL20S	ICONSL15M	ICONSL20M	ICONSL15F	ICONSL20F
King Single	1980x1050	ICONKS15S	ICONKS20S	ICONKS15M	ICONKS20M	ICONKS15F	ICONKS20F
King Single, Long	2030x1050	ICONKSL15S	ICONKSL20S	ICONKSL15M	ICONKSL20M	ICONKSL15F	ICONKSL20F
Double	1980x1350	ICOND15S	ICOND20S	ICOND15M	ICOND20M	ICOND15F	ICOND20F
Double, Long	2030x1350	ICONDL15S	ICONDL20S	ICONDL15M	ICONDL20M	ICONDL15F	ICONDL20F
Queen	1980x1530	ICONQ15S	ICONQ20S	ICONQ15M	ICONQ20M	ICONQ15F	ICONQ20F
Queen, Long	2030x1530	ICONQL15S	ICONQL20S	ICONQL15M	ICONQL20M	ICONQL15F	ICONQL20F
Split Queen	1980x760	ICONSQ15S	ICONSQ20S	ICONSQ15M	ICONSQ20M	ICONSQ15F	ICONSQ20F
Split Queen, Long	2030x760	ICONSQL15S	ICONSQL20S	ICONSQL15M	ICONSQL20M	ICONSQL15F	ICONSQL20F
King	1980x1850	ICONK15S	ICONK20S	ICONK15M	ICONK20M	ICONK15F	ICONK20F
King, Long	2030x1850	ICONKL15S	ICONKL20S	ICONKL15M	ICONKL20M	ICONKL15F	ICONKL20F

Manufactured in Australia, Fully Customisable

The Icon, fully manufactured by Forté Healthcare in Armidale, Australia, is available in custom sizes with adjustable length, width, and thickness, using premium Australian and European materials.





The Mattress Makes The Bed™

Address

10 Ribbon Gum Ave, Armidale NSW 2350

Phone

1300 731 163

Email

sales@fortehealthcare.com.au

Website

fortehealthcare.com.au

Socials



© Copyright 2026 Forté Healthcare

